

Overhead suspended shield & LED Lamp

350/45/LED130F-001

(40 cm x 50 cm 2.0 mm Pb equivalent glass shield & LED 130F examination lamp)

Effective protection for the upper body


This 40 cm x 50 cm shield uses 2.0 mm lead equivalent glass for effective x-ray protection to the eyes and thyroid. The robust yoke assembly provides stability and protection to the window which carries a 0.5 mm lead equivalent flexible curtain beneath to help further reduce scatter radiation. The LED 130F lamp is a compact yet practical examination light with a light output of 70,000 lux.

Features: LED lamp

- > Light intensity at 1 metre: 70,000 Lux
- > Colour rendering index Ra 95
- > Focusable light field size: 14 cm to 25 cm
- > Electronic dimming
- > Working distance: 70 cm to 140 cm
- > Colour temperature: 4300 Kelvin
- > Temperature increase at head area: 0.5°C
- > Total power consumption: 28 W
- > Light source: 19 LED's
- > Life span of over 40,000 hours
- > Diameter of lamp head: 33 cm
- > Detachable sterilizable handle included

Features: Shield


- > 40 cm x 50 cm x-ray shield
- > 2.0 mm lead equivalent glass
- > 0.5 mm lead equivalent flexible curtain

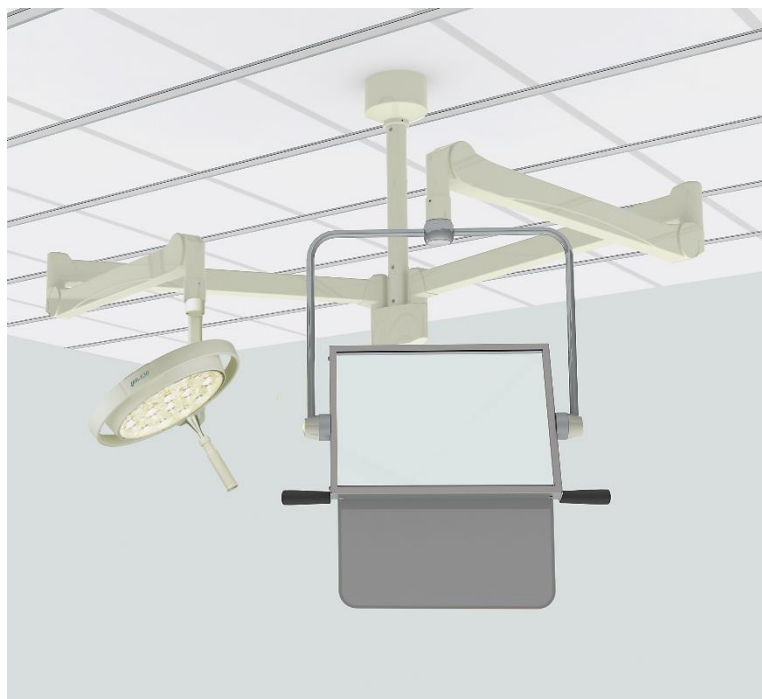

Both the shield and lamp are supported on two-part suspension arms, each attached to an 80 cm long ceiling column (other lengths available – see over). The column can be attached to a load bearing, or suspended ceiling. Different installation aids are available to suit differing site conditions (available separately – see over).

Alternatively, the ceiling column can be directly attached to the movable carriage of a 2.5 m or 4.2 m ceiling track (available separately) with the lamp cables contained and managed within a separately supplied flexible energy chain. These options enable the complete assembly to move into and out of the working area as required.

Technical Specification

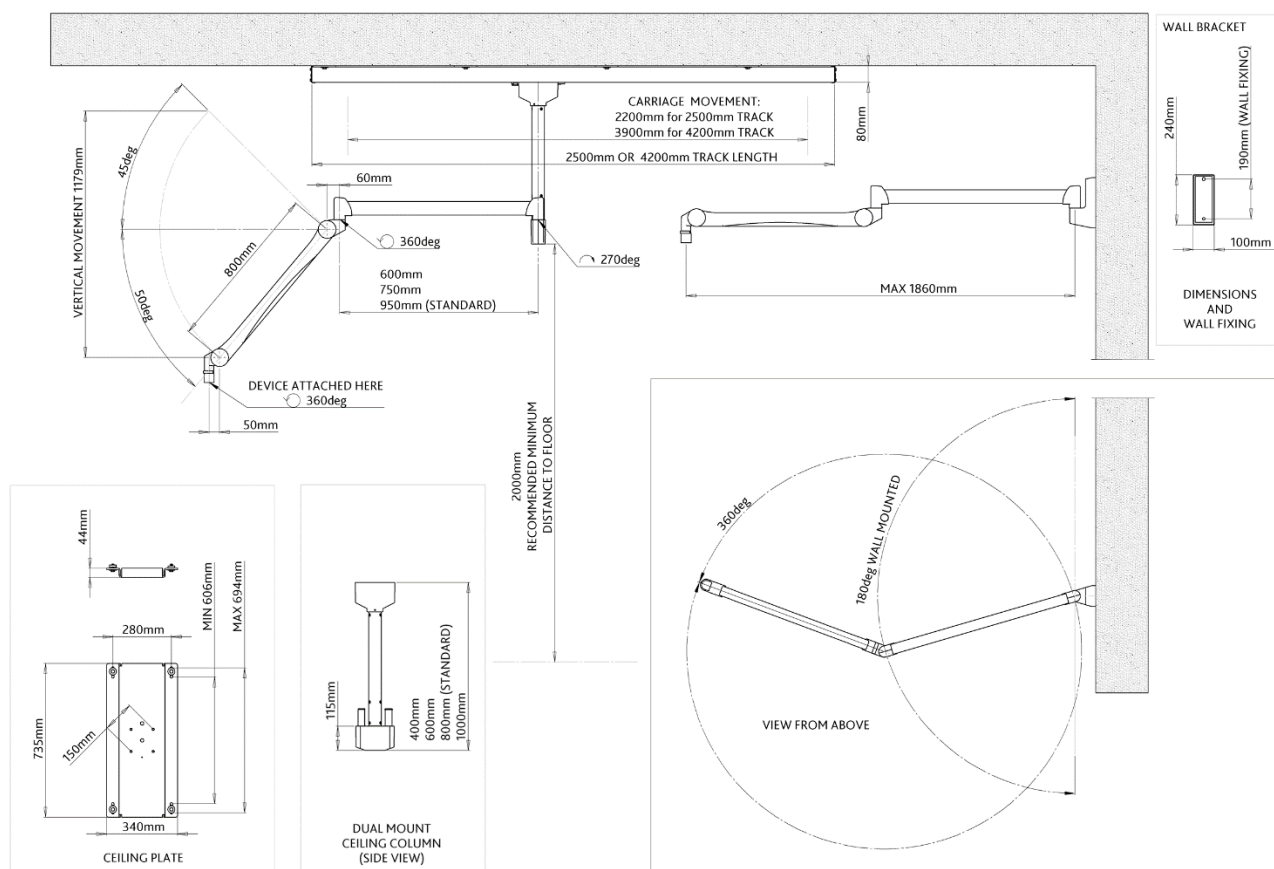

Equipment weights & part numbers

Included:

Power supply in housing
LED 130F lamp head: **2 kg** (approximately)
40 cm x 50 cm lead glass shield: **12 kg** (approximately)
2 x Suspension arm: **9 kg** (each)
80 cm dual mount ceiling column: **12.8 kg**

Optional equipment:

Dual mount ceiling column:

40 cm (#300/COLUMN/40): **10.2 kg**
60 cm (#300/COLUMN/60): **11.5 kg**
100 cm (#300/COLUMN/100): **14.1 kg**

Ceiling track & carriage:

2.5 m (#3001/TC/250): **50 kg**

4.2 m (#3001/TC/420): **72 kg**

Track cable management:

2.5 m (#3001/TC/250/EC)

4.2 m (#3001/TC/420/EC)

Additional track carriage:

For use in 2.5 m or 4.2 m track (#3001/3.25): **4.6 kg**

Optional installation aids:

- Ceiling column anchorage plate #300/P/004 (shown above attaching to ceiling rails with separately supplied ceiling column)
- Ceiling void installation kits #300/P/012 & #300/P/014

